

Prepared for:

**Centers for Medicare & Medicaid Services and Office of the
National Coordinator for Health Information Technology**

**CMS Alliance to Modernize Healthcare
Federally Funded Research and Development Center**

**Technical Authority for the Unified Clinical Quality
Improvement Framework**

Bonnie User Guide

Version 6.0

November 29th, 2020

Record of Changes

Version	Date	Author / Owner	Description of Change
0.1	April 15, 2014	Andre Quina / MITRE	Initial draft
0.1.1	July 18, 2014	Cristen Jones / MITRE	Update for Bonnie v 1.1
0.2	April 13, 2015	Cristen Jones / MITRE	Update for Bonnie v 1.3
1.4	January 13, 2016	Pace Ricciardelli / MITRE	Update for Bonnie v 1.4
1.5	November 18, 2016	David Connolly and Lizzie DeYoung / MITRE	Update for Bonnie v 1.5
1.6	January 4, 2017	Kristian Mulcahy and Chris Tohline / MITRE	Update for Bonnie v 1.6
1.6.5	June 8, 2017	Chris Tohline / MITRE	Update for Bonnie v 1.6.5
2.0	November 14, 2017	Rose Almonte, Lizzie Charbonneau, Chris Tohline, Jenn Reed / MITRE	Update in preparation for Bonnie v 2.0
2.1	May 3, 2018	Lizzie Charbonneau / MITRE	Update for Bonnie v 2.1
2.1.2	May 31, 2018	Lizzie Charbonneau / MITRE	Update for Bonnie v 2.1.2
2.2	July 13, 2018	Cole Springate-Combs / MITRE	Update for Bonnie v 2.2
2.3	August 23, 2018	James Bradley / MITRE	Update for Bonnie v 2.3
3.0	October 16, 2018	James Bradley / MITRE	Update for Bonnie v 3.0
3.1	November 15, 2018	James Bradley / MITRE	Update for Bonnie v 3.1
3.2	May 21, 2019	James Bradley / MITRE	Update for Bonnie v 3.2
4.0	August 27, 2019	James Bradley / MITRE	Update for Bonnie v 4.0
4.2	October 19th, 2020	Nicole Hunter / SemanticBits	Update for Bonnie v 4.2
6.0-Beta	November 29th, 2020	Nicole Hunter / SemanticBits	Updates for Bonnie FHIR v 6.0-Beta

Table of Contents

1. Introduction	1
1.1 Background	1
1.2 Purpose	1
1.3 Application Description	1
2. User Account Creation	3
2.1 Login Page	3
2.2 Creating a New User	3
2.3 Resetting a Password	4
2.4 Account Management	5
3. Measure Dashboard	6
3.1 Overview	6
3.2 Loading a New Measure	7
3.2.1 Value Set Options	9
3.3 Updating an eCQM	10
3.4 Creating Synthetic Test Records	11
4. Measure Results View	12
4.1 Overview	12
4.2 Measure Logic	15
4.3 Creating a New Test Record	16
4.4 Calculation Results	16
4.5 Editing a Test Record	18
4.6 Cloning a Test Record	19
4.7 Deleting a Test Record	19
4.8 Updating a Measure	19
4.9 Deleting a Measure	19
4.10 Updating Measurement Period	19
5. Building a Patient Test Record	21
5.1 Overview	21
5.2 Building a Synthetic Patient	23
5.3 Building the Patient History	24
5.4 Incremental Calculation	28
6. Feedback and Support	29
7. Frequently Asked Questions	31
Acronyms	33

List of Figures

Figure 1. Bonnie Login Page	5
Figure 2. Account Registration Page	6
Figure 3. Password Reset Page	6
Figure 4. Account Management Page	7
Figure 5. Measure Dashboard View	8
Figure 6. New Measure Dialog	9
Figure 7. Value Set Options – Profile	10
Figure 8. Value Set Options - Release	11
Figure 9. Updating Measure Dialog	12
Figure 10. Measure View	15
Figure 11. Measure Terminology and Overlapping Value Sets	16
Figure 12. Expanded Results View	18
Figure 13. Logic Calculation Highlight – Passing Results	18
Figure 14. Logic Calculation Highlight – Expanded Result	19
Figure 15. Logic Calculation Highlight – Failing Results	20
Figure 16. Change Measurement Period	21
Figure 17. Patient Builder View	22
Figure 18. Continuous Variable Measures Expected Populations	24
Figure 19. Building Patient History, Including Edit Clinical Element View	25

1. Introduction

1.1 Background

Bonnie is a software tool that allows electronic clinical quality measure (eCQM) developers to test and verify the behavior of their eCQM logic. The main goal of the Bonnie application is to reduce the number of defects in eCQMs by providing a robust and automated testing framework. The Bonnie application allows measure developers to independently load measures that they have constructed using the Measure Authoring Tool (MAT). Developers can then use the measure metadata to build a synthetic patient test deck for the measure from the clinical elements defined during the measure construction process. By using measure metadata as a basis for building synthetic patients, developers can quickly and efficiently create a test deck for a measure. The Bonnie application helps measure developers execute the measure logic against the constructed patient test deck and evaluate whether the logic aligns with the intent of the measure.

Bonnie has been designed to integrate with the nationally recognized data standards the Centers for Medicare & Medicaid Services (CMS) quality reporting programs use for expressing eCQM logic for machine-to-machine interoperability. This integration provides enormous value to the

eCQM program and federal policy leaders and stakeholders. The Bonnie tool verifies that the new and evolving standards for eCQMs used in the CMS quality reporting programs are flexible and can be implemented in software.

Bonnie was also designed to provide an intuitive and easy-to-use interface based on feedback from the broader measure developer community. A key goal of the Bonnie application is to deliver a user experience that provides an efficient and intuitive method for constructing synthetic patient records for testing and validating eCQMs.

Finally, the Bonnie software is freely available via an Apache 2.0 open source license. The CMS Electronic Health Records (EHR) Incentive Program makes all or parts of the Bonnie software available for inspection, verification, and even reuse by other government programs or federal contractors.

1.2 Purpose

The purpose of this document is to describe the functionality of the Bonnie web application that allows measure developers to test and verify the behavior of their eCQM logic. This document provides Bonnie users with step-by-step instructions for testing eCQMs by building synthetic patient records.

1.3 Application Description

The Bonnie application provides the capability to import measures defined following the Fast Healthcare Interoperability Resource (FHIR) Quality Measure Implementation Guide and Clinical Quality Language (CQL). The FHIR Quality Measure Implementation Guide specification provides the metadata and the CQL libraries provide the logic for calculating an eCQM. The Bonnie application can load an eCQM containing the FHIR and CQL describing a measure and programmatically convert both into an executable format that allows calculation of the measure directly from the specification.

The CMS MAT is the primary source for eCQMs used by the Bonnie application. Measure developers use the MAT to build eCQMs and export those measures as measure bundles containing FHIR and one or more CQL libraries. These measure bundles can be downloaded from the MAT and loaded into the Bonnie user interface for measure testing.

Once an eCQM has been loaded into the Bonnie application, a user can inspect the measure logic and then build synthetic test records and set expectations on how those test records will calculate against a measure. This capability to build synthetic test patient records, set expectations against those records, and calculate the measures using those patient records provides an automated and efficient testing framework for eCQMs. Through the Bonnie-supported eCQM testing framework, measure developers can more clearly understand the behavior of the measure logic and validate that the measure logic encodes their intent, as well as validate multiple iterations of measure updates against a test deck.

In addition, the development of a test deck as part of measure development provides benefits after the measures are finalized. For example, the test deck can demonstrate the intent of the measure through the patient examples in the deck. The systems in the test deck implement those measures with a means to validate their development through a base set of synthetic patient records with known expectations for calculating against the implemented measures.

2. User Account Creation

2.1 Login Page

Users require a valid account to access the Bonnie application. All eQMs loaded into the Bonnie application and all synthetic test patients are isolated by account. Therefore, users can only view, access, and modify data that they have loaded under their own account. Figure 1 shows the login screen for the Bonnie application. To log in, a user must provide the email address and password for a valid account.

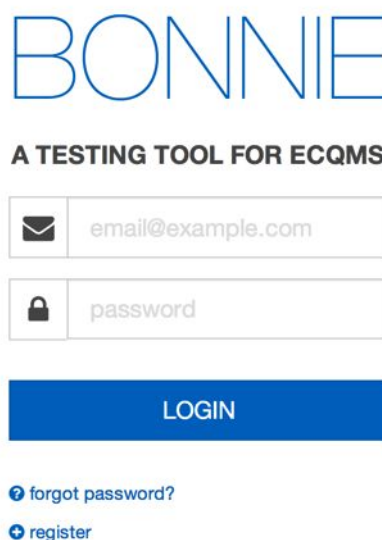
The image shows the Bonnie login page. At the top, the word "BONNIE" is displayed in a large, blue, sans-serif font. Below it, the text "A TESTING TOOL FOR ECQMS" is written in a smaller, bold, black font. There are two input fields: the first is for an email address, with a placeholder "email@example.com" and an envelope icon; the second is for a password, with a placeholder "password" and a lock icon. Below these fields is a blue button with the text "LOGIN" in white. At the bottom, there are two links: "forgot password?" with a question mark icon and "register" with a plus icon.







Figure 1. Bonnie Login Page

2.2 Creating a New User

A user can create a new account by clicking the “register” link on the login page. The register link brings the user to the account creation page shown in Figure 2. A user can create a new account by filling out the fields in the registration form and clicking the “REGISTER” button. Once an account has been created, your account will then go to approval. You will receive an email from the bonnie help desk asking for some additional information and then the account can be approved. After the users account has been approved the user can log in to the Bonnie application with the email address and password specified as part of account creation.

BONNIE

REGISTER

 first name	 last name
 email@example.com	 (000) 000-0000
 password*	 confirm password*

Must be at least 8 characters long and contain characters from at least two of: lowercase letters, uppercase letters, numbers, and special characters.

Terms and conditions for use of Bonnie.

Performance measures and related data specifications (the "Measures") are copyrighted by the noted quality measure providers as indicated in the applicable Measure. Coding vocabularies are owned by their copyright owners. By using the Measures, a user ("User") agrees to these Terms of Use. Measures are not clinical guidelines and do not establish a standard of medical care and quality measure providers are not responsible for any use of or reliance on the Measures.

THE MEASURES AND SPECIFICATIONS ARE PROVIDED "AS IS" WITHOUT ANY WARRANTY OF ANY KIND, AND ANY AND ALL IMPLIED WARRANTIES ARE HEREBY DISCLAIMED, INCLUDING ANY WARRANTY OF NON-INFRINGEMENT, ACCURACY, MERCHANTABILITY AND FITNESS FOR A PARTICULAR PURPOSE. NO QUALITY MEASURE PROVIDER, NOR ANY OF THEIR TRUSTEES, DIRECTORS, MEMBERS, AFFILIATES, OFFICERS, EMPLOYEES, SUCCESSORS AND/OR ASSIGNS WILL BE LIABLE TO YOU FOR ANY DIRECT, INDIRECT, SPECIAL, EXEMPLARY, INCIDENTAL, PUNITIVE, AND/OR CONSEQUENTIAL DAMAGE OF ANY KIND, IN CONTRACT, TORT OR OTHERWISE, IN CONNECTION WITH THE USE OF THE BONNIE TOOL OR THE MEASURES, EVEN IF ADVISED OF THE POSSIBILITY OF SUCH LOSS OR DAMAGE.

The Measures are licensed to a User for the limited purpose of using, reproducing and distributing the Measures with the Bonnie Tool as delivered to User, without any modification, reproduction and distribution of the Measures is prohibited. The user agrees to the Measures available on the Bonnie Tool and Terms of Use of the Measures.

☐ I agree to the above Terms and Conditions of use*

[forgot password?](#) [sign in](#)

CANCEL **REGISTER**


Figure 2. Account Registration Page

2.3 Resetting a Password

If a user forgets a password or an account is locked, the user can reset the password using the password reset page shown in Figure 3. This page is accessed from the “forgot password?” link on the login page (Figure 1). On the password reset page, the user provides the email address associated with the account and then presses the “Send” button. This action sends an email to the registered email address for the account to facilitate resetting the user’s account password.

EVERYONE FORGETS.

Just enter your email address, and we'll send you a link to reset your password.

	email@example.com
---	-------------------

SEND


[+ register](#)

[+ sign in](#)

Figure 3. Password Reset Page

2.4 Account Management

After logging into the application, the user can change the information associated with the account by accessing the account management page shown in Figure 4. The account management page can be opened by clicking the “Account” link in the application header shown at the top of Figure 4. From this page, users can change the information provided during the registration process and select a new password for their account.

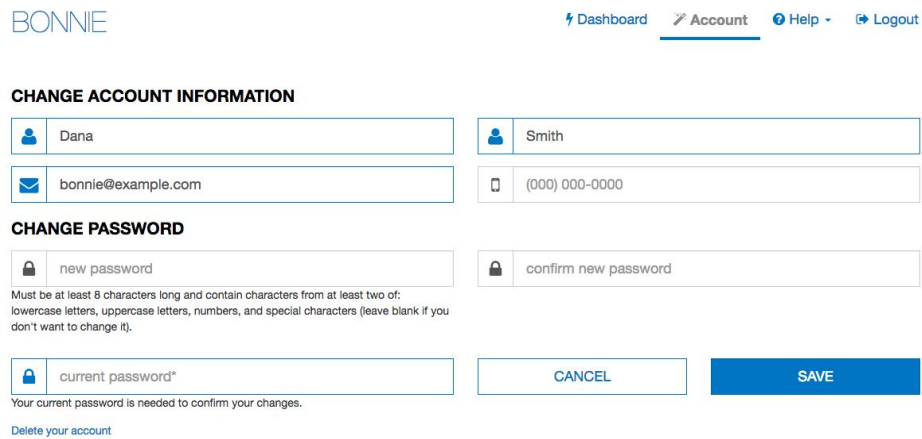

The screenshot shows the 'BONNIE' application header with navigation links: Dashboard, Account (active), Help, and Logout. The main content area is titled 'CHANGE ACCOUNT INFORMATION' and contains two columns of input fields. The left column has fields for 'First Name' (Dana), 'Email' (bonnie@example.com), 'New Password' (new password), and 'Current Password' (current password*). The right column has fields for 'Last Name' (Smith), 'Phone' ((000) 000-0000), and 'Confirm New Password' (confirm new password). Below the password fields, a note states: 'Must be at least 8 characters long and contain characters from at least two of: lowercase letters, uppercase letters, numbers, and special characters (leave blank if you don't want to change it).'. At the bottom, there are 'CANCEL' and 'SAVE' buttons, and a link to 'Delete your account'.

BONNIE [Dashboard](#) [Account](#) [Help](#) [Logout](#)

CHANGE ACCOUNT INFORMATION

CHANGE PASSWORD

Must be at least 8 characters long and contain characters from at least two of:
lowercase letters, uppercase letters, numbers, and special characters (leave blank if you
don't want to change it).

Your current password is needed to confirm your changes.

[Delete your account](#)

Figure 4. Account Management Page

3. Measure Dashboard

3.1 Overview

The Measure Dashboard page, as shown in Figure 5, is the initial page presented to users when they log into the application. The Measure Dashboard shows the calculation status of each eCQM loaded into the system. The calculation status shows how many patient records have been built for the measure, whether the measure is currently passing or failing, and how many patients are passing or failing for each measure.

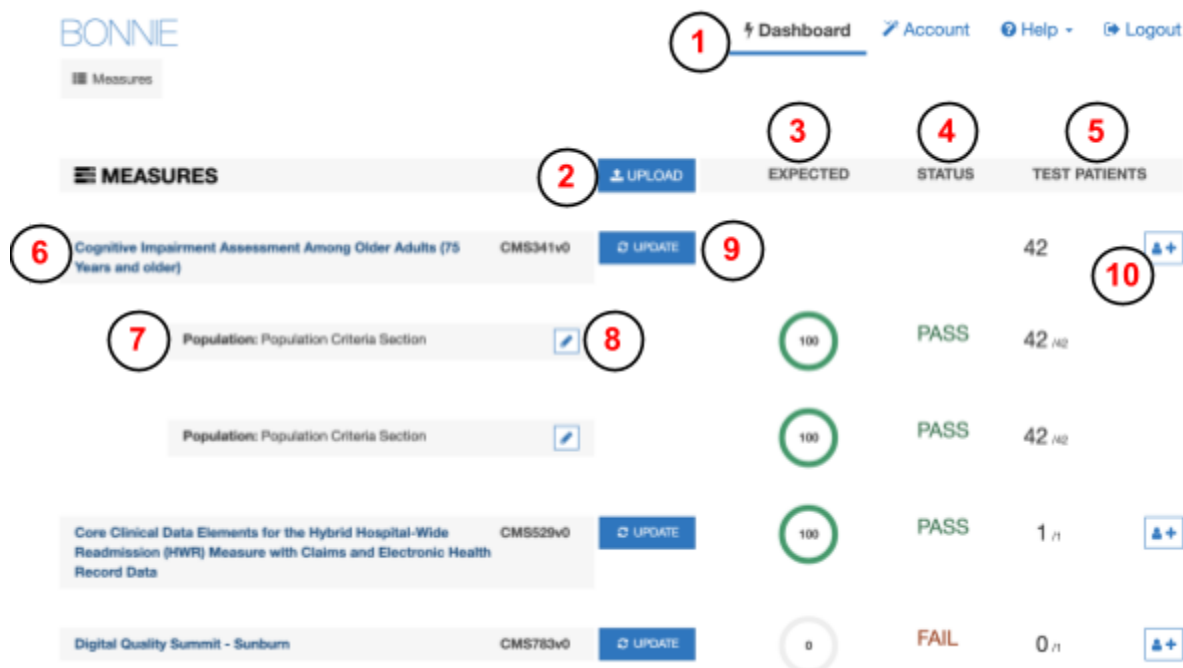


Figure 5. Measure Dashboard View

The Measure Dashboard View employs the following user interface (UI) elements (indicated by their item numbers in Figure 5):

1. Header – Allows the user to navigate to different parts of the application, access account information, send a support email (Contact), and log out of the application.
2. Upload Button – Allows the user to upload a new eCQM.
3. Expected Column – Displays the percentage of passing patients for the eCQM.
4. Status Column – Displays the status of the eCQM (New, Pass, Fail).
5. Test Patient Column – Displays the number of patients passing out of the total number of patients.
6. Measure Title – Displays the title for the eCQM and allows navigation to the Measure Dashboard View.

7. Subpopulation and stratification titles – Displays the titles for subpopulations or stratifications of an eCQM.
8. Edit Title Button – Allows the user to rename a subpopulation or stratification.
9. Update Button – Allows the user to update a previously loaded eCQM.
10. Add Patient Button – Allows the user to start building a new patient for an eCQM.

3.2 Loading a New Measure

When users log into the system for the first time, no eCQMs will be associated with the account. The user's first step is to load an eCQM into the account to begin testing the eCQM with the Bonnie application. The New Measure dialog, as shown in Figure 6, prompts the user to upload an eCQM.

The user may either upload an eCQM now or close this display and upload an eCQM later. After one or more eCQMs have been uploaded, this page displays the current set of eCQMs loaded into the system along with the subpopulations and stratifications associated with the eCQMs. The Measure Dashboard allows users to navigate to the details of individual eCQMs, upload a new eCQM, or update the definition of an existing eCQM. To upload a new eCQM, click the "Upload" button (item #3 in Figure 5) on the Measure Dashboard, which opens the New Measure Dialog shown in Figure 6.

The screenshot shows the 'New Measure' dialog box with the following components and numbered callouts:

- 1** MEASURE DATA: Includes a 'BROWSE' button and a file name 'CMS144_v5_4_Artifacts.zip'.
- 2** VSAC API KEY: Includes a text input field for 'API Key' and a link to 'Register for VSAC'.
- 3** VALUE SETS: Includes tabs for 'Profile' and 'Release', a dropdown menu for 'Latest eCQM {eCQM Update 2020-05-07}', a 'Draft' checkbox, and an 'Override options with measure defined versions' checkbox.
- 4** TYPE: Includes tabs for 'Eligible Professional' and 'Eligible Hospital'.
- 5** CALCULATION: Includes tabs for 'Patient Based' and 'Episode of Care'.
- 6** Includes a checkbox for 'Include Supplemental Data Element Calculations'.
- 7** CLOSE: A button to close the dialog.
- 8** LOAD: A button to load the measure.

Figure 6. New Measure Dialog

After clicking the “Upload” button, the user needs to fill out the information requested from the New Measure Dialog. These fields are outlined below and are indicated by their corresponding item numbers in Figure 6.

1. Browse – Opens a file browser enabling the user to select a MAT measure package to upload into Bonnie.
2. VSAC API Key – Field for entering the user’s VSAC API key to enable Bonnie to load the value sets included within the eCQM. You can find your API key by going to <https://uts.nlm.nih.gov/>, signing in to your UTS account, and clicking “My Profile.” If you do not already have an API key, you can generate one by clicking “Edit Profile” then checking the box that says “Generate new API Key” then clicking “Save Profile.”
3. Value Set Options – Options to define which value sets should be downloaded from VSAC. More information about these options is included in subsection 3.2.1.
4. Type – Option to define if the eCQM is an Eligible Professional/Eligible Clinician or Eligible Hospital/Critical Access Hospital eCQM.
5. Calculation – Option to define if the eCQM is a Patient-based eCQM or an Episode of Care-based eCQM.
6. Include Supplemental Data Element Calculations – Option to define if supplemental data elements should be included in calculations and displayed in the supplemental data element section in the eCQM logic. **Please note that in Beta SDE calculations will not function properly.**
7. Close – Closes the New Measure Dialog without uploading the eCQM.
8. Load – Uploads the specified eCQM using the selected options.

The action of clicking the “Load” button in the New Measure Dialog uploads the eCQM to the application for processing. Once eCQM loading is complete, the application directs the user to the Measure Dashboard (Figure 5) with the new measure available.

3.2.1 Value Set Options

Bonnie allows users to download value sets associated with an eCQM based on either a VSAC profile or VSAC release. Bonnie also allows users to download draft value sets (if “Profile” is selected) and/or override their selections with value set versions specified within the CQL logic itself as shown in Figure 7.

VALUE SETS

1 Profile Release

2 Latest eCQM <eCQM Update 2020-05-07>

3 ☐ Override options with measure defined versions

4 ☒ Draft

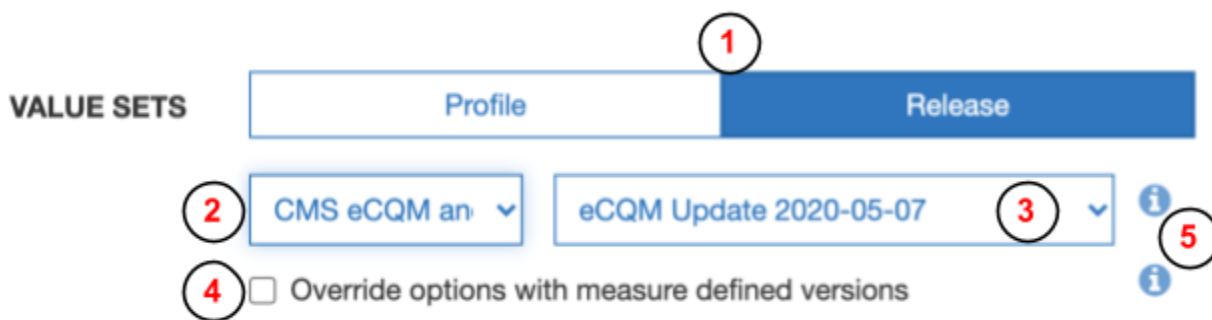
5 i

Figure 7. Value Set Options – Profile

Figure 7 shows the Value Set Options section (item #3 in Figure 6) with the default selections. The default selections are “Profile” with the “Latest eCQM” profile selected and “Draft” checked. These default selections have the same behavior as the default selections in Bonnie v2.0. The UI elements included when “Profile” is selected are as follows:

1. Profile – When Profile is selected, Bonnie queries the Value Set Authority Center (VSAC) for value sets calculated by an expansion profile. Expansion profiles are instructions for calculating value set content based on specific code system versions and retired legacy codes.
2. Profile dropdown selector – Allows the user to select an expansion profile to use for retrieving value sets. The default selection is “Latest eCQM <<{profile}>>” where {profile} is the VSAC profile referenced by “Latest eCQM”. In this screen shot, the referenced profile is “eCQM Update 2020-05-07”, so the display shows “Latest eCQM <<eCQM Update 2020-05-07>>”.
3. Override options with eCQM-defined versions – It is possible for an eCQM to contain value set versions within the CQL itself. This option allows the user to override the selected option to use the eCQM-defined value set versions. If this option is selected and a version for a value set is **not** defined within the measure, then Bonnie will use the selected profile to retrieve the appropriate value set version.
4. Draft – When this option is selected, Bonnie will retrieve the draft value set expansions for the selected expansion profile. Note that if the user selects draft and is not an assigned author or steward group member for one or more of the value sets included in the measure, the measure upload will fail. Draft is selected by default.
5. Information buttons – Contains additional information regarding the value set options.

Figure 8 shows the Value Set Options section (item #3 in Figure 6) with “Release” selected.

**Figure 8. Value Set Options - Release**

The UI elements included when “Release” is selected are as follows:

1. Release – When Release is selected, Bonnie queries VSAC for the published value sets associated with the specified release.

2. Program dropdown selector – Allows the user to select the program that includes the release they would like to use.
3. Release dropdown selector – Allows the user to select the appropriate release.
4. Override options with measure-defined versions – It is possible for an eCQM to contain value set versions within the CQL itself. This option allows the user to override the selected option to use the measure-defined value set versions. If this option is selected and a version for a value set is **not** defined within the eCQM, then Bonnie will use the selected release to retrieve the appropriate value set version.
5. Information buttons – Contains additional information regarding the value set options.

3.3 Updating an eCQM

Once the eCQM has been loaded, the testing process may identify issues with the eCQM. When issues are identified, the logic must be updated in the MAT to resolve these issues. Alternatively, the measure could be updated in the MAT as part of an Annual Update. After an eCQM has been updated in the MAT, it may be necessary to update that eCQM in Bonnie for testing. To update an eCQM, follow these steps:

1. Click the “Update” button (item #10 in Figure 5) on the Measure Dashboard, which displays the “Update Measure” dialog as shown in Figure 9.
2. Select a new MAT export zip with the updated eCQM definition.
3. Click the “Load” button to load the new version of the eCQM.

Follow-Up Care for Children Prescribed ADHD Medication (ADD) [✕]

MEASURE DATA	BROWSE CMS136v9.1.005__MATv5_8_Artifacts.zip
VSAC API KEY	API Key ? Register for VSAC
VALUE SETS	<div>Profile Release</div> <div>Latest eCQM «eCQM Update 2020-05-07» ?</div> <div><input checked="" type="checkbox"/> Draft ?</div> <div><input type="checkbox"/> Override options with measure defined versions ?</div>
TYPE	
CALCULATION	Patient Based

CLOSE LOAD

Figure 9. Updating Measure Dialog

3.4 Creating Synthetic Test Records

Once a set of eCQMs has been loaded into the Bonnie application, users can start building test patients for the eCQMs. To build a test patient from the Measure Dashboard (Figure 5), click the

“Add Patient” button (item #12). This action opens the patient builder screen as shown in Figure 17. For more information, please refer to Section 5, *Building a Patient Test Record*.

After the user creates synthetic test patients for eQMs, the Measure Dashboard displays summary calculation results of the patients associated with each eQM loaded by the user (Figure 5). As shown in Figure 5, UI elements numbered 4, 5, and 6 on the Measure Dashboard provide the summary results for each eQM. In the “Expected” column (item #5), the UI displays the percentage of patients associated with the eQM whose calculated values meet the expectations set for the patient. The “Status” column (item #5) indicates the current state of the measure—whether it is passing, failing, or new. An eQM is in the passing state if all associated patients are meeting expectations. If one or more patients are not meeting expectations, then the eQM is in a failing state. Any eQM is in the new state when there are no patients associated with it. Finally, the “Test Patients” column (#6) displays the number of passing patients over the total number of patients as a fraction. The measure view, as represented in Figure 10 in Section 4, shows more detailed results for an eQM.

4. Measure Results View

4.1 Overview

As shown in Figure 10, the Measure View page displays the detailed information, associated patients, and calculation results for a single eCQM. In this view, users can add new patients to an eCQM, update an eCQM, and delete an eCQM. To access the Measure View, click on the “Measure Title” link (item #7) on the Measure Dashboard, as shown in Figure 5.

The Measure View page presents the title and description of the eCQM, the logic for the eCQM, and a measure complexity indicator. The page lists the current set of test patients associated with the eCQM in the column on the right-hand side of the page. This column displays the summary calculation results for the patient test deck associated with the eCQM along with the list of patients and the individual calculation results for each. This page also shows the results of calculation for a single patient along with an overlay of the calculation results on the eCQM logic.

The following UI elements appear at the top of the Measure View page (indicated by their item numbers in Figure 10):

1. Measure Title – Displays the title and description of the eCQM.
2. Measurement Period – Displays the current Measurement Period and allows the user to update the Measurement Period.
3. Measure Actions – Allows the user to delete or update an eCQM definition.
4. Measure Subpopulations or Stratifications – Allows access to different subpopulations or stratifications in the eCQM.
5. Measure Logic – Displays a representation of the logic for the eCQM.
6. Logic Highlighted With Coverage – Displays the logic for the eCQM, highlighting which lines of the eCQM are covered by the test patients.
7. Patient Actions – Allows the user to export patient records as QRDA or Excel.
8. Percent Successful – Displays the percent of patients currently meeting expectations for the eCQM.
9. Patients Passing Count – Displays the current number of patients meeting expectations over the total number of patients in the test deck for the eCQM.
10. Test Coverage – Displays the percentage of the measure logic that has been evaluated to true for the patient test deck. This measure provides a method for determining how much of the logic has been tested.
11. Add Patient – Allows the addition of a new patient to the test deck for this eCQM.
12. Failing Patient – Displays an example of a patient that is not currently meeting expectations for the eCQM.
13. Passing Patient – Displays an example of a patient that is meeting expectations.

14. Patient Name – Displays the name given to the patient.
15. Patient Status – Displays PASS or FAIL to indicate if the patient is meeting expectations.
16. Expand Patient Results Button – Displays the calculation details of a patient. This display will show the expected and actual values for the patient against the eCQM.
17. Show Result – If a patient has been selected with the Expand Patient Results Button, this displays the result of the logic evaluated on the selected patient.
18. Show All Results – If a patient has been selected with the Expand Patient Results Button, this displays the result of each eCQM logic block evaluated on the selected patient.

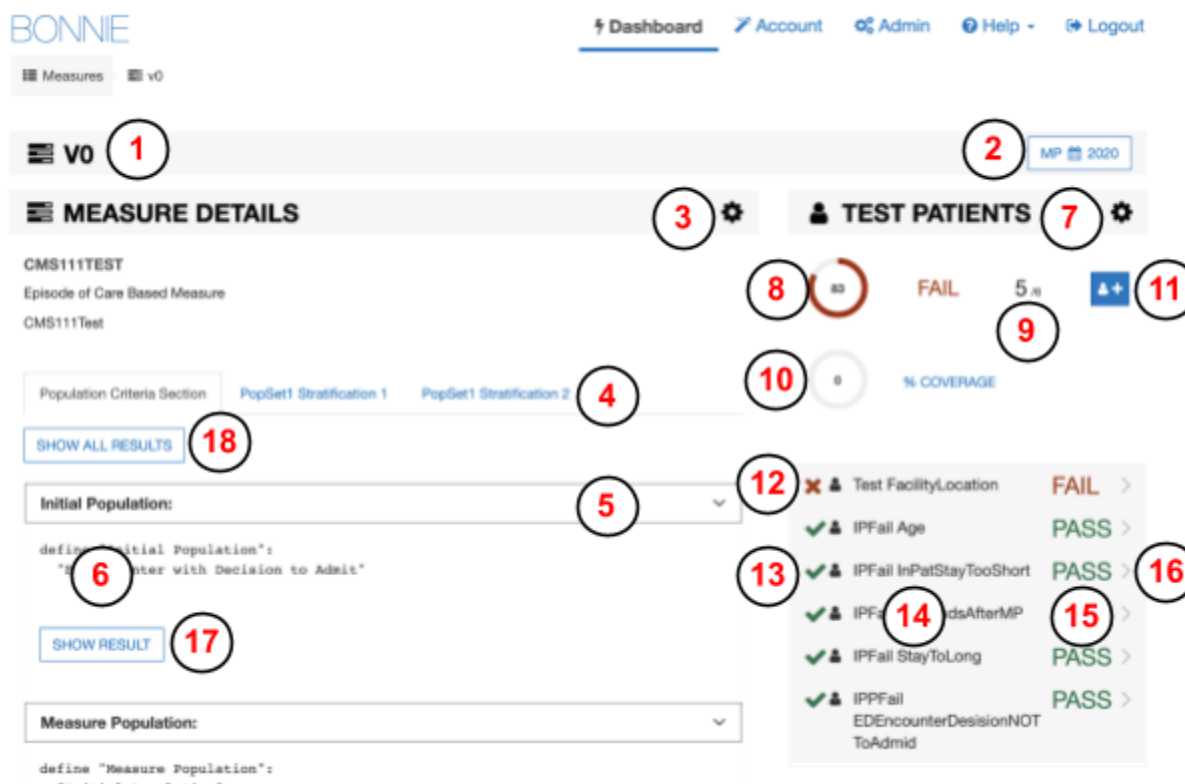


Figure 10. Measure View

The following UI elements (indicated by their item numbers in Figure 11) appear at the bottom of the Measure View Page:

1. Terminology – Displays the value sets from the eCQM and referenced libraries. Clicking on any one of these value sets will show the list of codes associated with that value set.
2. Overlapping Value Sets – Displays value sets in the eCQM and referenced libraries that have shared codes.

1 TERMINOLOGY		
Antipsychotic Medications	2.16.840.1.113883.3.464.1003.196.12.1255	Draft
Bipolar Disorder	2.16.840.1.113883.3.67.1.101.1.128	Draft
Birth date	Direct Reference Code	N/A
Dead	Direct Reference Code	N/A
Emergency Department Visit	2.16.840.1.113883.3.117.1.7.1.292	Draft
Encounter Inpatient	2.16.840.1.113883.3.666.5.3001	Draft
Ethnicity	2.16.840.1.114222.4.11.837	Draft
Huntington's Disease	2.16.840.1.113883.3.464.1003.105.12.1032	Draft
ICU Admission or Transfer	2.16.840.1.113883.17.4077.3.2040	Draft
Inpatient Psychiatric NDNQI	2.16.840.1.113762.1.4.1098.3	Draft
ONC Administrative Sex	2.16.840.1.113762.1.4.1	Draft
Payer	2.16.840.1.114222.4.11.3591	Draft
Race	2.16.840.1.114222.4.11.836	Draft
Schizophrenia	2.16.840.1.113883.3.464.1003.105.12.1104	Draft
Threat to themselves or others	2.16.840.1.113883.3.464.1003.195.12.1020	Draft
Tourette's Syndrome	2.16.840.1.113883.3.464.1003.105.12.1030	Draft
2 OVERLAPPING VALUE SETS		
Encounter Inpatient 2.16.840.1.113883.3.666.5.3001	Inpatient Psychiatric NDNQI 2.16.840.1.113762.1.4.1098.3	1 code
Inpatient Psychiatric NDNQI 2.16.840.1.113762.1.4.1098.3	Encounter Inpatient 2.16.840.1.113883.3.666.5.3001	1 code

Figure 11. Measure Terminology and Overlapping Value Sets

4.2 Measure Logic

The left-hand side of the Measure View contains the measure title, description, and a representation of the logic. The representation of the eCQM logic is similar to the human-readable display for the eCQM provided in the MAT measure exports.

The eCQM logic section can be used to verify that the eCQM logic was properly loaded from the measure bundle without the loss of any logical conditions. The logic can also be used to evaluate the nature of the calculation of a test patient against the logic (please refer to subsection 4.4, *Calculation Results*) and to visualize the test coverage of the eCQM logic.

4.3 Creating a New Test Record

To create a new test record, begin by clicking the “Add Patient” button (item #11) in the Measure View (Figure 10). This action opens the Patient Builder (shown in Figure 17). After creating the test patient record, the application returns the user to the Measure View where the user can evaluate the results of calculating the patient against the eCQM.

4.4 Calculation Results

The user can calculate the test patient against the logic of the eCQM in the Measure View (Figure 10). High-level results are calculated automatically when the Measure View is loaded. These high-level results appear under the test patients section on the right-hand side of the Measure View. These results include the percent of patients passing (item #8), the test coverage (item #10), individual statuses for each patient (items #12-13), and the passing patient count (item #9).

In addition to these high-level results, detailed results can be displayed for an individual patient by clicking the “Expand Patient Results” button (#16). This displays a table of expected and actual results for the patient covering each population listed in the measure and shows whether that population currently meets or fails expectations.

The Expanded Results View employs the following UI elements (as indicated by their item numbers in Figure 12):

1. Failing Population – A population for which the patient fails.
2. Passing Population – A population for which the patient passes.
3. Population Column – A list of the population types.
4. Expected Value – The user-defined expected value for the population.
5. Actual Value – The calculated value for that population.
6. Edit Patient Button – Allows editing of the selected patient.
7. Clone Patient Button – Allows cloning of the selected patient.
8. Delete Patient Button – Allows deleting of the selected patient.

IPPass

AssessmentPerformed

FAIL

	<div>3</div> Population	<div>4</div> Expected	<div>5</div> Actual
<div>1</div>	IPP		<input type="checkbox"/>
<div>2</div>	DENOM	<input type="checkbox"/>	<input type="checkbox"/>
	DENEX	<input type="checkbox"/>	<input type="checkbox"/>
	NUMER	<input type="checkbox"/>	<input type="checkbox"/>


6

7

8

Figure 12. Expanded Results View

Clicking the “Expand Patient Results” button also displays the patient results calculated against each line of logic. This information is displayed in the measure logic section of the view by highlighting each logic clause. As shown in Figure 13, Figure 14, and Figure 15, a green highlight (accompanied by a solid underline) indicates a passing result for the logic calculation, while a red highlight (accompanied by a double-underline) indicates a failing result over the applicable lines of text.

Definitions:


```

define "Inpatient Encounters":
  ["Encounter, Performed": "Encounter Inpatient"] Encounter
  where Global."LengthInDays"(Encounter.relevantPeriod)<= 120
  and Encounter.relevantPeriod ends during "Measurement Period"
  
```

SHOW RESULT

Figure 13. Logic Calculation Highlight – Passing Results

Figure 13 shows the results of a single patient calculated against the eCQM logic. The highlighting of the eCQM logic text for calculation is intended to provide a clearer understanding of why a patient is calculating to a specific population, and uses the following indications of status:

- Logic evaluated to TRUE – Green highlighting along with a solid underline
- Logic evaluated to FALSE – Red highlighting along with a double underline

The results of the calculation in Figure 13 are that the patient aligns with the logic of the initial patient population (IPP). The highlighting of the logic in Figure 13 indicates that every AND condition evaluated to true and at least one condition from each OR evaluated to true. Based on this calculation, the “Inpatient Encounters” definition evaluates to true for the patient.

The “Show Result” button at the bottom of Figure 13 can be used to expand the result, as shown in Figure 14. To collapse the results, click the “Hide Result” button.

The screenshot shows a user interface for eCQM logic calculation. At the top, a green header bar contains a checkmark icon, the text "Initial Population:", and a dropdown arrow. Below this, the logic definition is displayed with green highlighting and solid underlines: `define "Initial Population":` and `"Encounter With Age Range and Without VTE Diagnosis or Obstetrical Conditions"`. A blue button labeled "HIDE RESULT" is positioned to the left of a large grey box. The grey box contains a JSON-like structure representing the patient's encounter data:

```
[{
  account: null,
  appointment: null,
  basedOn: null,
  class: null,
  classHistory: null,
  contained: null,
  diagnosis: null,
  episodeOfCare: null,
  encounter: null,
  ...
}]
```

Figure 14. Logic Calculation Highlight – Expanded Result

In Figure 14, the bottom grey box displays information regarding the result of the eCQM logic definition evaluated on the selected patient. In the example shown, the “Inpatient Encounters” definition looks for encounters in a certain time range. The grey box will display all encounters that meet that criteria; in this example, the selected patient has one encounter that meets the criteria.

In Figure 15, the population is highlighted in red, indicating that the logical statements evaluate to false. Some clauses in this figure are green, indicating they are true; however, the false female

statement causes the AND logic to be false, and therefore, the results calculated in Figure 15 indicate that the patient is neither included in the IPP nor the Denominator.

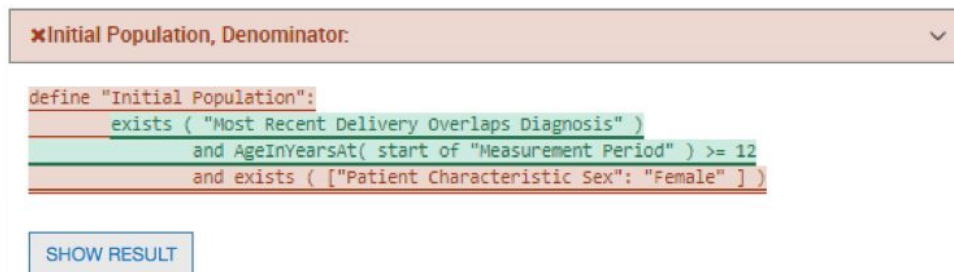


Figure 15. Logic Calculation Highlight – Failing Results

4.5 Editing a Test Record

The user can edit a test patient from the Measure View (Figure 10) by clicking the “Edit” button. The user accesses the “Edit” button for a patient (item #6 in Figure 12), accessible after expanding the patient result. Clicking the “Edit” button opens the Patient Builder (as shown in Figure 17) with the data populated for that patient. Once a patient record has been edited and saved, the application returns the user to the Measure View.

4.6 Cloning a Test Record

The user can clone a test patient from the Measure View (Figure 10) by clicking the “Clone” button (item #7 in Figure 12) to the immediate right of the “Edit” button, accessible after expanding the patient result. This action opens the Patient Builder (Figure 17) with the data populated for the patient being cloned. The difference between editing and cloning a patient is that the clone process creates a new patient based on an existing patient, while the edit process updates the data for an existing patient. Once a patient record has been cloned, edited, and saved, the application returns the user to the Measure View.

4.7 Deleting a Test Record

The user can delete a test patient from the Measure View (Figure 10) by clicking the “Delete” icon (item #8 in Figure 12) to the immediate right of the “Clone” button, accessible after expanding the patient result. Deleting a patient requires a two-step process for confirmation. Once a user deletes a patient record, the action cannot be undone. To delete a patient record, the user initially clicks the “Delete” button. A second “Delete” button is then displayed. The user must click the second “Delete” button to confirm the deletion of the patient.

4.8 Updating a Measure

The user can update a measure by clicking the “Update Measure” button, which is accessed by clicking the “Measure Actions” icon (item #3 in Figure 10). The “Update Measures” button displays the Update Measure Dialog (shown in Figure 9), which allows the user to specify a new zip package for an eCQM exported from the MAT. Once the new measure package has been

entered, the user clicks the “Load” button, which updates the eCQM definition and returns the user to the Measure View with the updated measure definition.

4.9 Deleting a Measure

The user can delete an eCQM from the Measure View (Figure 10) by clicking the “Delete” icon for an eCQM. To access the “Delete” icon, the user clicks the “Measure Actions” icon (item #3 in Figure 10). A user cannot undo the deletion of an eCQM. To delete an eCQM, the user initially clicks the “Delete” icon. A second “Delete” icon is then displayed. The user must click the second “Delete” icon to confirm the deletion of the eCQM. Deleted eCQMs are no longer displayed on the Measure Dashboard.

4.10 Updating Measurement Period

The user can modify the Measurement Period associated with the Measure from the Measure View (Figure 10) by clicking the “Measurement Period” button (item #2). Once the user clicks the button, a dialog will appear (Figure 16) which allows the user to input a new Measurement Period. Currently Bonnie does not allow you to shift all dates on the patient history by the same number of years that the measurement period will shift. This checkbox is disabled currently.

Change Measurement Period ×

	CURRENT		NEW
YEAR	<input type="text" value="2020"/>	→	<input type="text" value="2020"/>
PATIENTS	<input type="checkbox"/> Shift dates in patient data. <i>This option shifts all dates on all patients in this measure by the number of years being changed. Feb 29th in leap years will be changed to Feb 28th in non-leap years.</i>		

Figure 16. Change Measurement Period

5. Building a Patient Test Record

5.1 Overview

The Patient Builder view, as shown in Figure 17, allows the addition and editing of clinical data for a synthetic test patient record. The user accesses the Patient Builder view by clicking the “Add Patient” button (item #10) on the Measure Dashboard (Figure 5), or by clicking the “Add Patient,” “Edit,” or “Clone” buttons from the Measure View (Figure 10).

The Patient Builder View employs the following UI elements (as indicated by their item numbers in Figure 17):

1. Patient Name – Allows the entry of a first and last name for the patient record.
2. Patient Characteristics – Allows the definition of characteristics data for the patient.
3. Measure Information – Shows the description for the patient’s associated eCQM.
4. Expectations – Allows users to set the calculation expectation for each population of the eCQM.
5. Actions – Allows users to save or cancel a patient record.

Figure 17. Patient Builder View

The Patient Builder view also provides fields to either add new data or edit existing data for a patient. The patient’s first and last names can be defined using the associated edit fields.

Characteristics such as gender, birthdate, race, and ethnicity can be set in the characteristics section. Fast Healthcare Interoperability Resource elements can be added to the patient history section by dragging and dropping an individual element onto the patient history section.

In addition to defining the patient data, the Patient Builder view allows the user to set expectations on the patient using the “Expectations” section. Expectations represent how the user expects the patient to calculate against the eCQM. The “Expectations” section will be different based on the type of eCQM used to build the patient record. The “Expectations” section for patient-based eCQMs provides a check box for users to select the appropriate expected populations (e.g., numerator, denominator) within which the patient will be included. The “Expectations” section for episode of care-based measures allows users to select the number of episodes of care that are expected to be included in each population using a number picker. Continuous variable measures allow users to define the expected value(s) the measure is expected to calculate for the patient.

The Patient Builder’s logic section displays the logic of the eCQM against which the patient test record is constructed. As data is entered, the application continuously calculates the patient against the eCQM logic, and the results are displayed against the eCQM using the logic highlighting. Subsection 4.4 provides additional information about the descriptions of the logic highlighting technique based on calculation results shown in Figure 13 and Figure 15.

5.2 Building a Synthetic Patient

The first step in building a synthetic patient record is defining the patient name. The patient’s first and last name can be set in the “Patient Name” section (item #1 in Figure 17). The user can then add patient notes and define patient characteristics in the “Patient Characteristics” section (item #2). Characteristics for the patient include data such as birthdate, race, ethnicity, gender and living status.

After defining the patient’s name, the user defines the expectations for how the patient will behave and be calculated against the eCQM. The user sets the expectations for the patient in the “Expected Value” section (item #4). If the user expects the patient to align with the initial patient population logic, then the user should set that expectation here, in the “Expected Value” section (item #4). Expectations are used throughout the Bonnie application to determine if a patient is passing or failing against the eCQM.

If the expected results for the patient align with the actual results from calculating the patient against the eCQM logic, the patient passes. If the expected results do not match the actual calculated results, the patient fails. This system provides the capability to build patients and set expectations for those patients based on the intent of the eCQM. Once the eCQM is calculated, the passing state indicates that the eCQM logic aligns with the intent, while the failing state indicates that either the eCQM logic does not align with the intent or the synthetic patient was constructed improperly.

Note that in Figure 17, the depicted measure has a discrete expected value. Accordingly, the NUMER, DENOM, and DENEX options are shown. For continuous variable measures that rely on an aggregate calculation (e.g., a median time between ED arrival and ED departure), the OBSERV, MSRPOPL, and (optionally) MSRPOPLEX options are shown, as in Figure 18.

EXPECTED VALUE

IPP ☒ MSRPOPL ☒

OBSERV % MINS 50

Figure 18. Continuous Variable Measures Expected Populations

Although discrete measures offer only a checkbox for including patients in the NUMER population, continuous variable measures offer the ability to set a certain percentage or minute value for a given patient calculation. This value is only available if the patient is in the MSRPOPL population.

5.3 Building the Patient History

Figure 19 depicts two events in the patient's history—an encounter and a clinical summary—based on data elements from the eCQM. When a FHIR resource is added to the patient history, it becomes an event in the patient's history. When an event is first created, any default timing elements are added and it is associated with a default code path code from each value set associated with the source FHIR resource. These defaults as well as additional data can be edited by expanding the event. Figure 19 shows an example of an expanded event from the patient history that can be edited. By clicking the expand/collapse details icon (item #3), the user can edit the details of the element.

The following UI elements are shown in Figure 19:

1. Elements Section – Contains FHIR elements that the user can add to the patient history.
2. FHIR Element – Shows a condensed summary of a FHIR element that has been added to the patient history.
3. Expand/Collapse Details – Allows hiding or expanding the details of an element.
4. Default Timing Section – This section allows the user to add default timing elements. If a timing element is a choice type (it can be different timing types) it needs to be added in the attributes section. These fields are defaulted to today's date and the measurement period year.
5. Codes Section – Allows adding the primary code path to the element. Note that Bonnie will automatically find an appropriate code and add it.
6. Attributes Section – Allows adding fields to the element (i.e., Result, Diagnoses).
7. Extension Section – Allows the user to enter FHIR extension attributes. URL is required, value is optional. Value can be many different types. First select your type then fill in the value
8. Delete Button – Allows deleting an element from the patient history.
9. CQL Calculation Results – Shows the results for each of the populations.

10. Show Result – Displays the result of the logic evaluated on the patient.
11. Show All Results – Displays the result for each eCQM logic block evaluated on the patient.

The screenshot displays the 'Edit Clinical Element View' interface. On the left, a sidebar titled 'ELEMENTS' contains a list of categories: CARE PROVISION, CLINICAL SUMMARY, DIAGNOSTICS, FINANCIAL SUPPORT, MANAGEMENT, and MEDICATIONS. The main area is titled 'Management: Encounter: Encounter Inpatient'. It features a 'PERIOD - period' section with 'START' and 'END' date/time pickers. Below this is a 'CODES' section with a text input for 'SNOMEDCT: 183452005' and a 'Code System' dropdown. The 'ATTRIBUTES' section includes a 'status: "finished"' dropdown. The 'EXTENSIONS' section has a 'URL' and a 'VALUE' dropdown. At the bottom, there is a 'CLINICAL SUMMARY' section with a 'DATE' dropdown. On the right, a 'V0' section shows 'CQL Calculation Results' with a table of results. Below this are 'SHOW ALL RESULTS' and 'SHOW RESULT' buttons. The interface is annotated with numbered callouts: 1 points to the 'ELEMENTS' sidebar; 2 points to the 'CLINICAL SUMMARY' section; 3 points to the 'V0' section; 4 points to the 'PERIOD - period' section; 5 points to the 'CODES' section; 6 points to the 'ATTRIBUTES' section; 7 points to the 'EXTENSIONS' section; 8 points to the 'VALUE' dropdown; 9 points to the 'CQL Calculation Results' section; 10 points to the 'SHOW RESULT' button; and 11 points to the 'SHOW ALL RESULTS' button.

Figure 19. Building Patient History, Including Edit Clinical Element View

After the user defines the patient characteristics, the user builds a patient history from FHIR elements extracted from the eCQM. The available elements from the eCQM are organized by category and listed in the “Elements” section (item #1). The user can click a category to expand the list of available elements. Expanding this list allows the user to click and drag an individual element onto the patient history (item #2).

Several fields in the Edit Clinical Element View can be edited for an event in the patient history using the controls shown in Figure 19. These fields include the default timings of the event (items #4), codes (item #5), various attributes (item #6), and extensions (item #7). The timing elements can be set for an event by either typing into the text fields directly or by using the date/time pickers that are displayed when the field is selected. An undefined start time and end time can be set for the event (used for active or ongoing events) by selecting the “Undefined” checkbox. Selecting the “Undefined” check box clears the date/times, indicating that the event has not started/ended.

Attributes can be set for an event. Attributes such as ordinal, severity, discharge date/time, etc. can be added using the “Attributes” section. Attributes are added by selecting the type of the field (scalar, coded, or time) and entering a scalar value, selecting a value set, or entering a date.

Finally, events can be removed from the patient history by clicking the “Delete” icon (item #10). Deletion requires a two-step process. After initially pressing the “Delete” icon, the user is prompted to depress the “Delete” icon a second time to confirm the deletion.

5.4 Incremental Calculation

The final section of the Patient Builder View is the logic section. This section displays a representation of the logic for the eCQM against which the patient is constructed. The logic can be a reference to help describe the details of the data that should be added in building the synthetic patient.

The logic section continuously displays the results of calculating the patient against the eCQM by means of the logic highlighting described in Figure 13 and Figure 15. Any modification made to a patient triggers a recalculation of the patient against the eCQM, which updates the results of the calculation displayed by the logic highlighting. Therefore, while the user constructs the patient record, the user can inspect the behavior of the logic relative to the patient.

When the user has completed constructing a synthetic patient, the user clicks the “Save” button in the Patient Builder View. This action adds the patient to the test deck for the measure and returns the user to the Measure View. After the user creates the first patient, additional patients can be created from scratch or cloned from existing patients to extend the coverage of the test deck against the eCQM.

6. Feedback and Support

An issue tracker and feedback email list are available to support the resolution of issues and to answer questions related to the Bonnie application. The Bonnie issue tracker is available on the ONC Jira system at: <http://jira.oncprojectracking.org/browse/BONNIE>

When users encounter bugs in the Bonnie application, they should use the issue tracker to report bugs, ask questions, or to request new features. To add an issue, users must create a login account in the Jira system. Once an issue has been entered, the Bonnie team will review and prioritize it.

7. Frequently Asked Questions

Does Bonnie replace Cypress¹ or is it an alternative to Cypress for certification?

No. Bonnie is a testing tool for measure developers to test eCQMs as they are being authored, while Cypress is the CMS quality reporting program certification tool. Bonnie cannot be used for CMS quality reporting program certification for vendors.

Bonnie provides insight into the details of how eCQM logic behaves. This makes it useful to vendors, measure developers, and others in the eCQM space who would like to better understand the behavior of eCQMs. The Bonnie testing tool uses the same measure calculation engine as the Cypress tool. Therefore, Bonnie can be used to better understand the behavior of the eCQMs when implementing an eCQM prior to certification.

Can I export patient records from Bonnie?

Not at this time. Bonnie will support this in a future iteration.

Can I load patient records into Bonnie?

Currently, Bonnie does not support loading patient records into the tool. If you would like to calculate eCQMs using existing patients, the popHealth tool may be a better solution.

Does Bonnie automatically generate patient records?

Currently, Bonnie does not offer the capability to automatically generate patient records based on the eCQM logic. The goal of Bonnie is to allow the construction of carefully considered patient scenarios. These synthetic patient records are intended to have an expected result based on the intent of the eCQM. Automatically generating patient records would generate patients that align with the eCQM logic as written, but if there were issues in the eCQM logic, these patients might not align with the intent of the eCQM.

My patient does not match the logic of the Initial Population. Why is the patient passing?

A test in Bonnie is based on assigning an expected outcome to a synthetic patient record based on the intent of the eCQM. A patient can be constructed with the expectation that the patient calculated is to be included in the Initial Patient Population. If the patient actually calculates to the Initial Patient Population, then the test passes. If the patient does not calculate to the Initial Patient Population, then the test fails.

There must be an exact match between the patient's alignment and the expectation's needs. For instance, if the expectation is set for only the Initial Patient Population, and the patient calculates to the Initial Patient Population and the Denominator, the test fails. Similarly, a patient can be constructed with the expectation that the patient does not align with any of the eCQM populations. In this case, if the patient does not calculate into the Initial Patient Population, the test will pass. If a patient is constructed with the expectation set of not aligning with any of the eCQM populations and that patient calculates into the Initial Patient Population, that test will fail.

¹ Cypress is the rigorous and repeatable testing tool of Electronic Health Records (EHR) and EHR modules used in calculating electronic clinical quality measures used in CMS's quality reporting programs.

Where can I get help with Bonnie?

The Bonnie application has a Help menu in the header. Once you log into the application, you can use the help menu to send an email to the Bonnie feedback list to (a) access the Bonnie issue tracker to enter a bug, ask a question, or request a feature; and (b) access the user guide that provides step-by-step instructions on using the Bonnie tool.

What measure formats can I load into the Bonnie tool?

The Bonnie application can load CQL 1.4 measure packages from the Measure Authoring Tool using the FHIR R4 model. Note that when loading, you will need a National Library of Medicine (NLM) Value Set Authority Center (VSAC) account (<https://uts.nlm.nih.gov/license.html>) to download the value sets associated with the measure.

Do I have to be a measure developer to use the Bonnie tool?

No. Anyone can sign up for a Bonnie account using the “register” link on the login page.

Do I need to be a Measure Authoring Tool user to use the Bonnie tool?

Normally, no, but Bonnie requires you import a Measure Authoring Tool (MAT) bundle to be successful. Since no FHIR measures bundles have been posted to the CMS website yet you will need to be a current MAT user or receive an export from a current MAT on FHIR user.

Can Bonnie be used to calculate the results for a large number of patient records?

Bonnie is not designed to handle calculations for more than a few hundred patient records per measure. If you are interested in calculating eCQMs against a larger number of patient records, the popHealth tool may be a better solution.

Acronyms

Acronym	Definition
CMS	Centers for Medicare & Medicaid Services
CQL	Clinical Quality Language
CQM	Clinical Quality Measure
eCQM	Electronic Clinical Quality Measure
ED	Emergency Department
EH	Eligible Hospital
EP	Eligible Professional
FHIR	Fast Healthcare Interoperability Resources
HHS	Department of Health and Human Services
HTML	Hypertext Markup Language
IPP	Initial Payment Population
MAT	Measure Authoring Tool
NLM	National Library of Medicine
ONC	Office of National Coordinator for Health Information Technology
QRDA	Quality Reporting Document Architecture
UI	User Interface
UMLS	Unified Medical Language System
VSAC	Value Set Authority Center
XML	Extensible Markup Language